



AperiSpray: A Low-Cost 3D-Printed Paper Spray Ion Source Conversion Kit for Commercial Electrospray Housings

Lucas R. Abruzzi^{1,2}, Chris G. Gill^{1,2,3,4*}

¹Center for Health and Environmental Mass Spectrometry (CHEMS), Department of Chemistry, Vancouver Island University, Nanaimo, BC, Canada

²Department of Chemistry, University of Victoria, Victoria, BC, Canada

³Department of Environmental and Occupational Health Sciences, University of Washington, Seattle, WA, USA

⁴Canadian Institute for Substance Use Research, University of Victoria, Victoria, BC, Canada


Table of Contents

Table of Contents	2
I.0 Overview of Source	3
I.1 Schematic	3
I.2 Parts list	4
2.0 Source construction	5
2.1 Printing parts	5
2.1.1 Filament	5
2.1.2 Printer settings	5
2.1.3 Slicer settings	5
2.2 Source Disassembly	8
2.3 Soldering Electronics	13
2.4 Assembling the shaft assembly	16
2.5 Assembling the High-Voltage Termination.	21
2.6 Installing the strip mount	23
2.7 Installing assembly in the ESI housing	27
3.0 Opensource Strip Cartridge Construction	29
3.1 Cutting Paper Strips	29
3.2 Printing Strip Holders and Strip Holder Clips	29
3.3 Assembling Strip Cartridges	30
4.0 Operating Instructions	32
4.1 Positioning Calibration Procedure	32
4.2 Spray Solvent Delivery Set Up	33
4.3 Method Development	33
4.4 Sample Preparation	34
4.5 Loading and Analyzing a Strip Cartridge	34
4.6 Troubleshooting	36
4.7 Unloading Used Strip Cartridges	36

I.0 Overview of Source

The AperiSpray is a low-cost, simple to use assembly that converts existing Thermo Scientific IonMax and OptaMax NG heated electrospray ionization (H-ESI) housings into paper spray ion sources, making use of 3D-printed parts and off-the-shelf hardware. AperiSpray features instrument-controlled high-voltage and spray solvent delivery and can use either in-house-built or commercially available paper spray strip cartridges. Many laboratories that perform routine LC-MS analyses already have surplus H-ESI sources available for conversion, and those that don't can purchase inexpensive used or defective housings from third-party sellers, making AperiSpray a low-investment option for those wanting to experiment with PS-MS. Additionally, unlike the commercial PS-MS platform, AperiSpray is compatible with almost any Thermo Scientific mass spectrometer with an atmospheric pressure interface, allowing it to be coupled to high-resolution Orbitrap and time-of-flight instruments.

I.1 Schematic

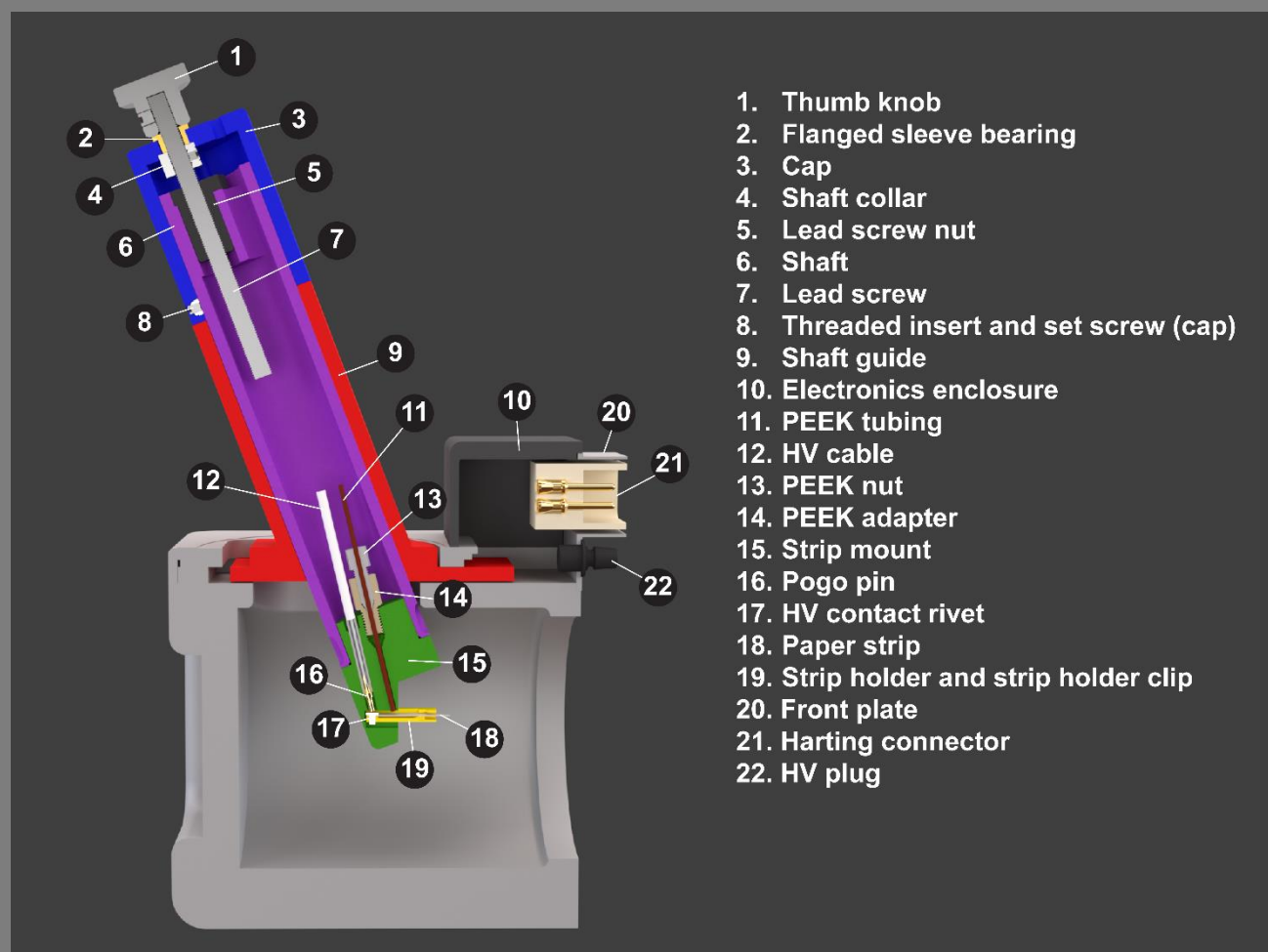

I.2 Parts list

Identifier	Qty	Source	Part No.	Cost /ea (CAD)	Material
Shaft guide	1	3D-printed	N/A	2.43	PLA-CF
Shaft	1	3D-printed	N/A	2.46	PLA-CF
Cap	1	3D-printed	N/A	1.46	PLA-CF
Electrical enclosure	1	3D-printed	N/A	1.60	PLA-CF
HV female connector	1	3D-printed	N/A	0.16	PLA-CF
Strip mount	1	3D-printed	N/A	0.72	PA6-GF
Strip holder + clip	N/A*	3D-printed	N/A	0.06	PA6-GF
Strip tray	1	3D-printed	N/A	0.73	PLA
Harting connector	1	Digikey.ca	1195-1808-ND	29.82	PC
2.15 k Ω resistor	1	Digikey.ca	56-SFR250000215IFR500CT-ND	0.20	Ceramic
22 M Ω resistor	1	Digikey.ca	BC4913CT-ND	0.74	Ceramic
20 AWG wire	4 cm	mcmaster.com	8054T184	0.02	Copper
Heat shrink	1	homedepot.ca	1000843481	25.41	Polyolefin
HV cable	70 cm	mcmaster.com	8296K17	25.67	Silicone
Crimp butt connector	1	grainger.ca	GGE4FRH9	0.13	Vinyl
Pogo pin	1	Digikey.ca	54-7949-0-15-20-09-14-11-0-ND	4.76	Brass alloy
Pen tube	1	grainger.ca	TOYGSM11-2	0.47	PP
1/4"-PP tubing	6 cm	mcmaster.com	5392K12	0.20	PP
Threaded insert	1	mcmaster.com	94100A110	0.70	SS
Set screw (cap)	1	mcmaster.com	92605A098	0.15	SS
Thumb knob	1	mcmaster.com	2577K21	12.57	SS
Flanged sleeve bearing	1	mcmaster.com	7815K11	5.17	Bronze
Lead screw	1	mcmaster.com	6350K694	81.44	SS
Shaft collar	1	mcmaster.com	6462K12	7.17	SS
Lead screw nut	1	mcmaster.com	6350K683	31.63	POM
Mounting screws	2	mcmaster.com	92196A151	0.11	SS
Set screw (strip-mount)	1	mcmaster.com	92765A005	0.33	Steel alloy
HV contact rivets	N/A*	aliexpress.com	3 mm M2 GB109	0.03	Aluminum
Chromatography paper	N/A*	www.fishersci.ca	057146	0.02	Cellulose
PEEK nut	1	coleparmer.ca	0201343	11.32	PEEK
PEEK nut adapter	1	coleparmer.ca	0201420	40.81	PEEK

*Denotes consumables

2.0 Source construction

2.1 Printing parts

2.1.1 Filament

For the cap, shaft, shaft guide, electronics enclosure, and strip tray, most standard FDM filament materials are acceptable. We prefer carbon-fiber PLA because it is dimensionally stable, easy to print and has a desirable surface texture. For the strip mount and strip holders, we recommend a material that is resistant to solvent to prevent interference and instrument contamination. In normal operation of the source, the strip mount should not come in direct contact with spray solvent, so PLA or ABS can be used if necessary, however in our experience, solvent overflow will occasionally occur while calibrating solvent deposition rate and volume, so a solvent-resistant material is highly recommended. Strip holders come in contact with solvent during normal operation. Some good options for strip holders (and strip mounts) include natural polyoxymethylene (POM), and natural PA6 Nylon. We recommend natural glass fiber reinforced PA6 nylon as it is both chemically resistant and easy to print on consumer grade printers, whereas POM suffers from poor bed adhesion and severe warping.


2.1.2 Printer settings

Printer settings are highly dependent on the filament material and printer. We recommend using your printer's default print settings for the filament chosen and adjusting as necessary.

2.1.3 Slicer settings


Shaft guide:

Print the shaft guide in an upright position. Tree support is recommended for the overhang in the loading cutaway.

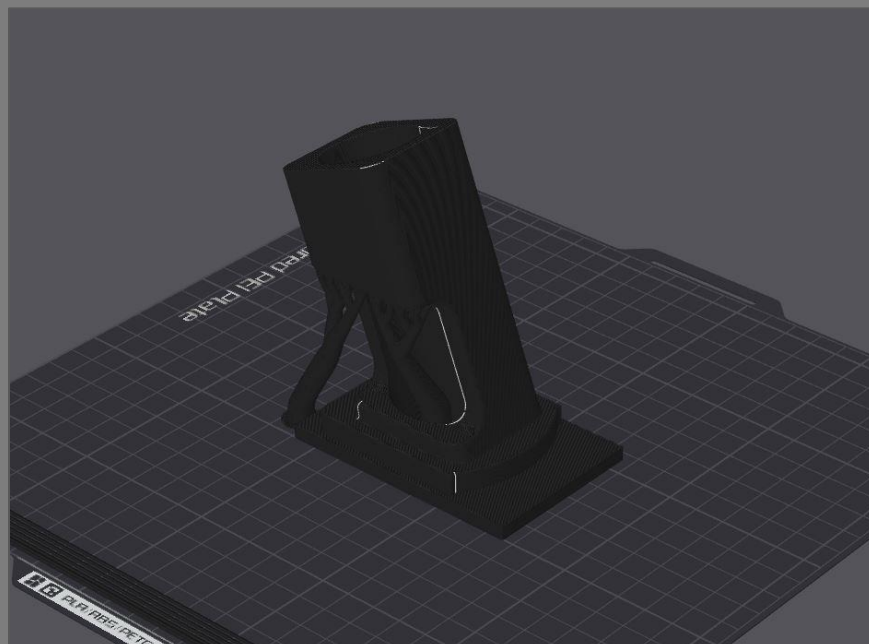

Shaft:

We recommend printing the shaft upright to ensure the surface that mates with the lead screw nut is flat. You will need support under the mounting bracket on the inside of the shaft.

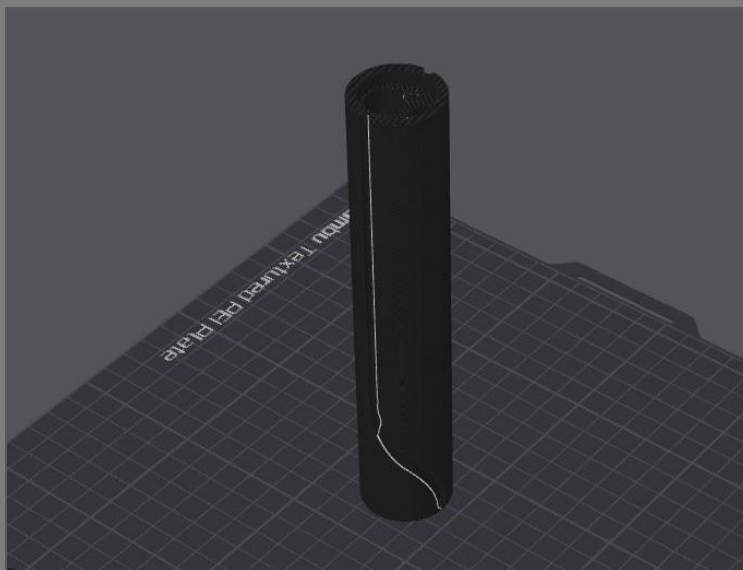


Cap:

The cap should be printed upside down with the top face against the build plate. No support is necessary. We recommend using a contrasting filament for the warning text and HV symbol if multi-filament printing is available.

Electronics enclosure:

Print the electronics enclosure with the backside down on the build plate. Support is only needed for the bottom lips.

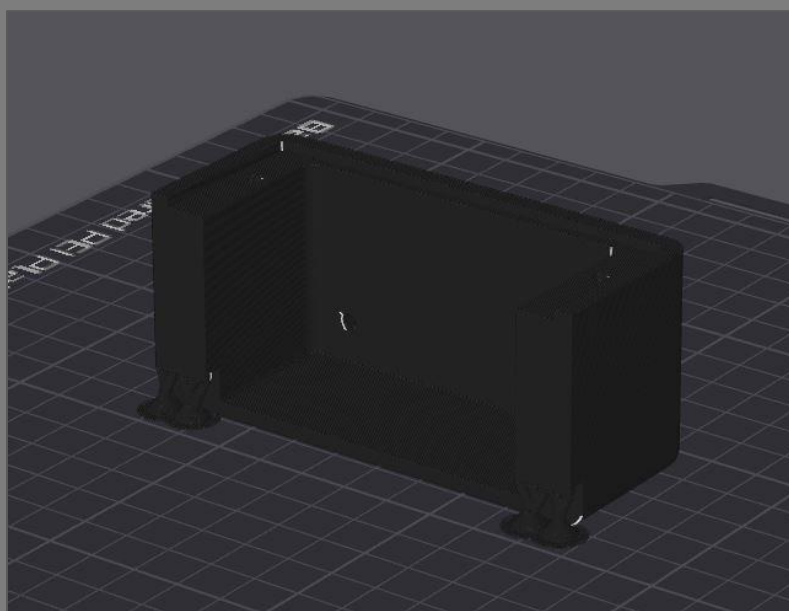

HV female connector:

Print with the threaded end up. No support is needed.

Strip mount (for open-source or commercial strip cartridges):

The strip mount is best printed upside down with support under the ridge, in the strip pocket, and in the threaded solvent hole.



Calibration Strips:


Calibration strips are for calibrating the strip position relative to the inlet of the mass spectrometer. They are inserted into the strip pocket on the strip mount and extend to where the opening of the inlet should be relative to the strip mount. By default, the calibration strips extend 4 mm further from the strip mount than actual paper strips would to create a 4 mm gap between paper strips and the inlet when used to calibrate positioning. If a larger or smaller gap is desired, the FreeCAD file for the calibration strips can be modified by increasing or decreasing the "Inlet_gap" variable in the variable set. Print calibration strips out of whatever filament you have on hand. They are best printed in an upright position with the backside against the print bed. Support is only necessary on the calibration strip for the commercial strip mount.

2.2 Source Disassembly


Remove the 4 screws holding the front plate on the source and pull the front plate away from the housing.




Remove the metal cover protecting the electronics connections.

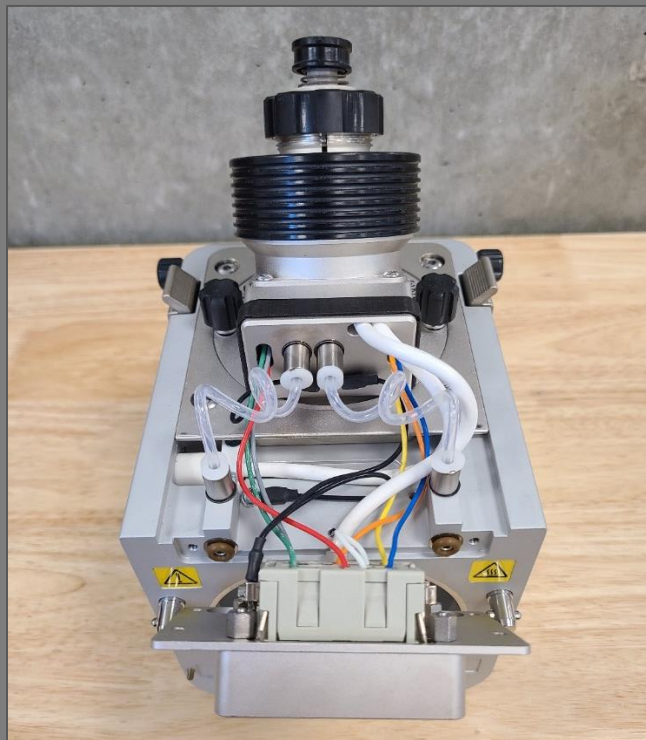


Disconnect gas lines by pushing down on quick connect tab while pulling up on gas tubing. Remove the retainer knobs for the API spray insert and heater assembly and the 4 socket head screws that mount the spray insert and heater assembly to the housing.

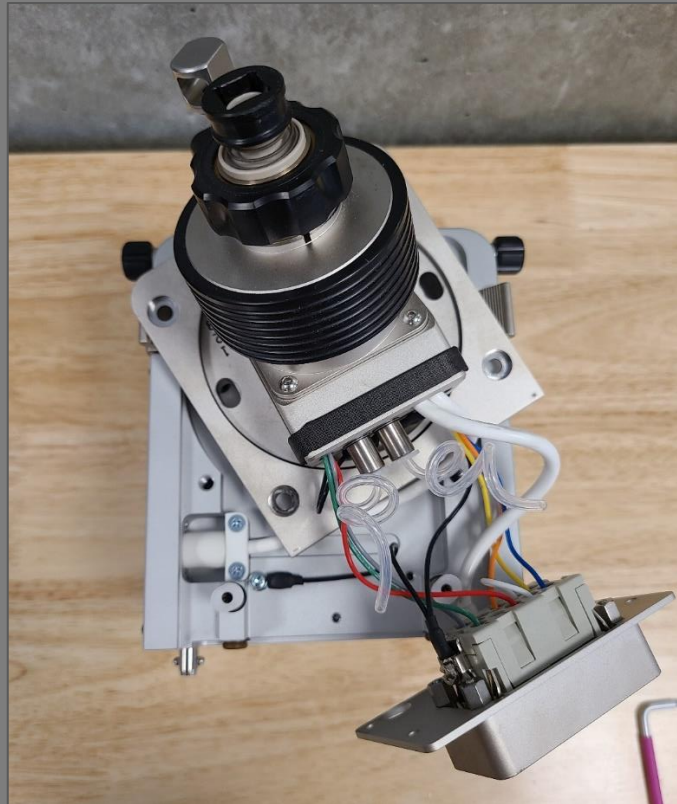


Remove the 3 socket head screws and the side panel on the right side of the housing.

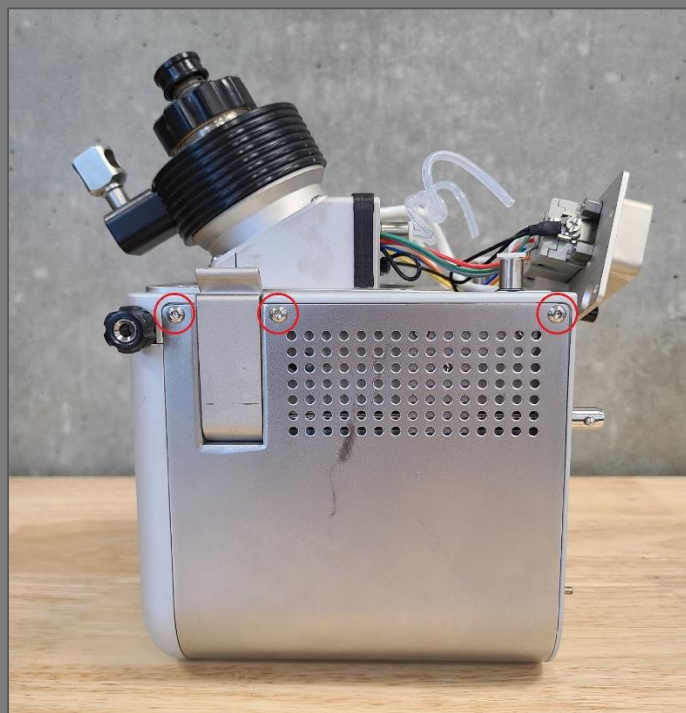


Remove the 4 socket head screws and metal plate covering the APCI corona discharge needle high-voltage connection.

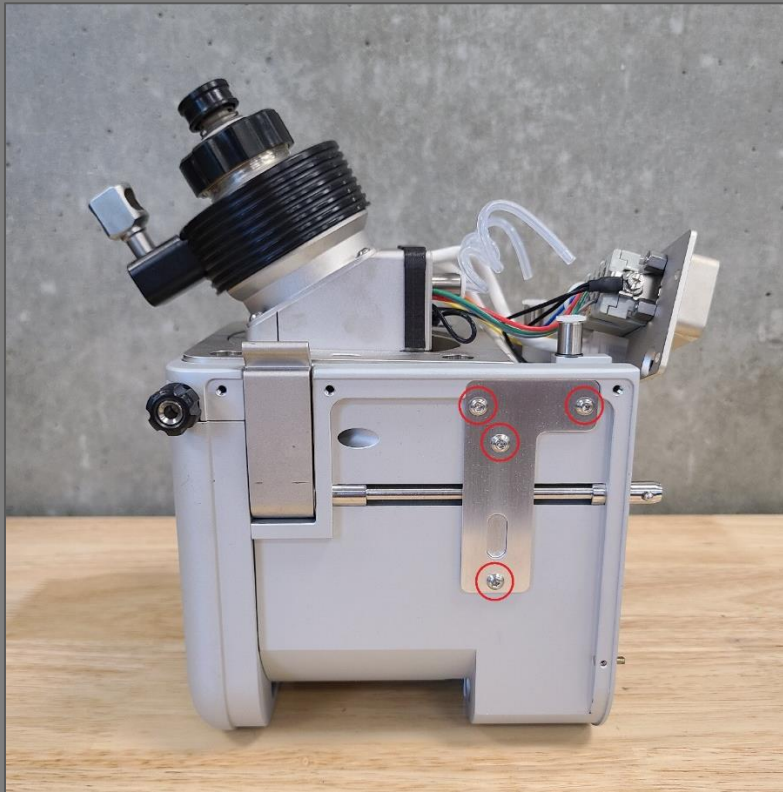


Remove the 2 screws and the plastic mount housing the high-voltage corona discharge needle connector and adjacent grounding connection. The right-side panel can be replaced.

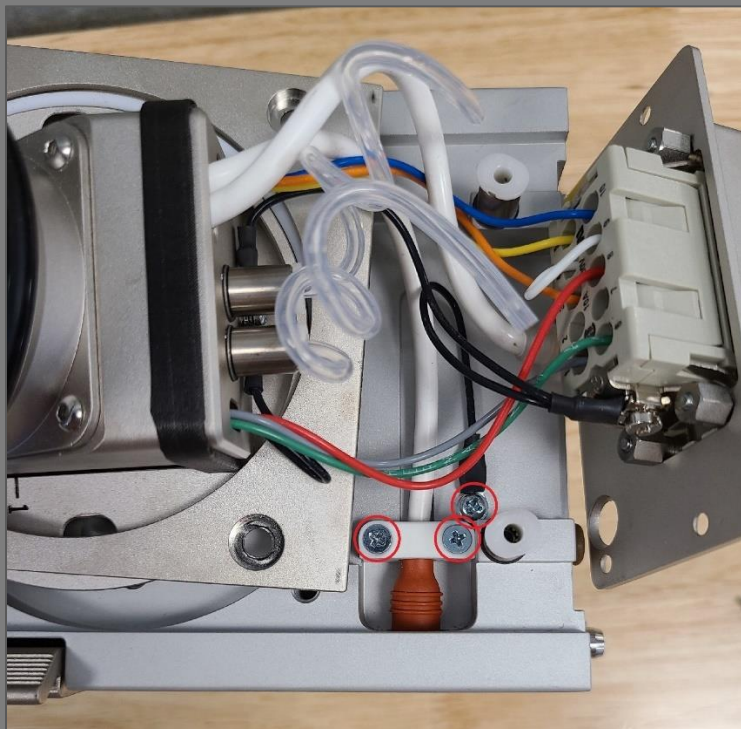


Remove the API spray insert and the 4 socket head screws and the heater assembly cover.

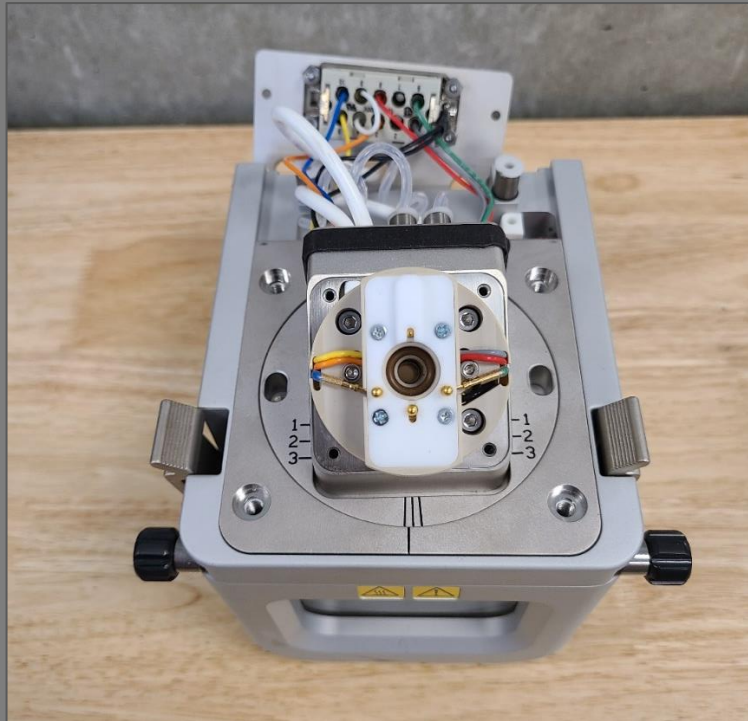


Remove 3 socket head screws holding on the heater assembly and carefully pull the heater assembly away from the mounting body while feeding wire.

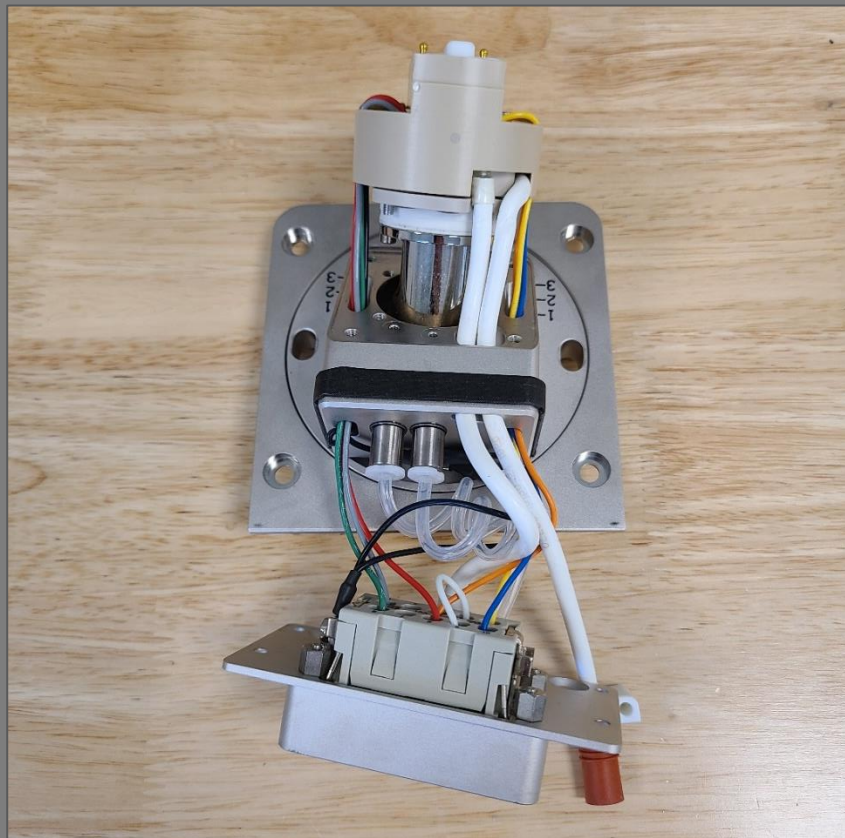

Carefully slide the positioning plates up and over the heater assembly and electrical connections to separate them from the electronic components

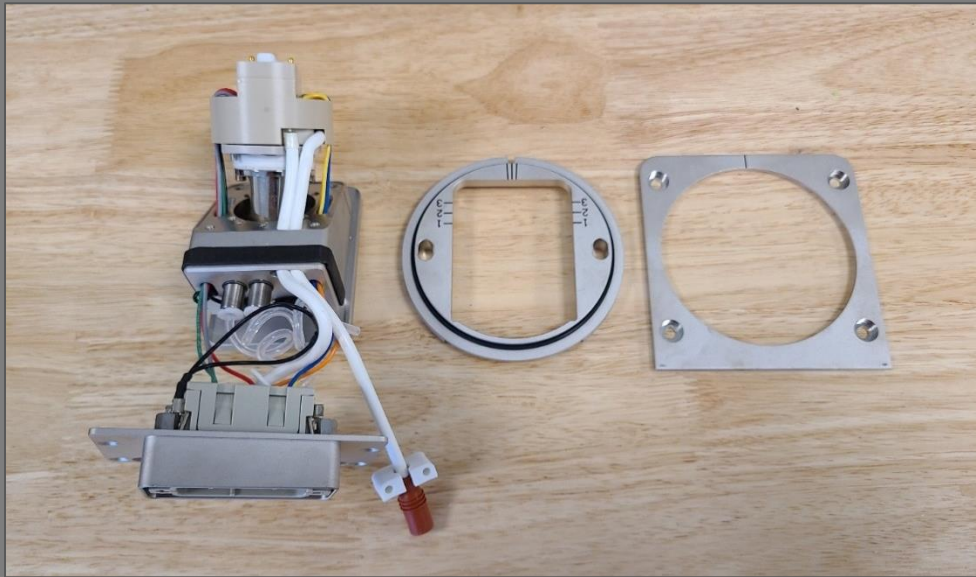


Loosen and remove the 4 nuts retaining the IO-pin Harting connector to the front plate and unscrew the transparent square nut connecting the HV wire to the HV plug on the front plate.

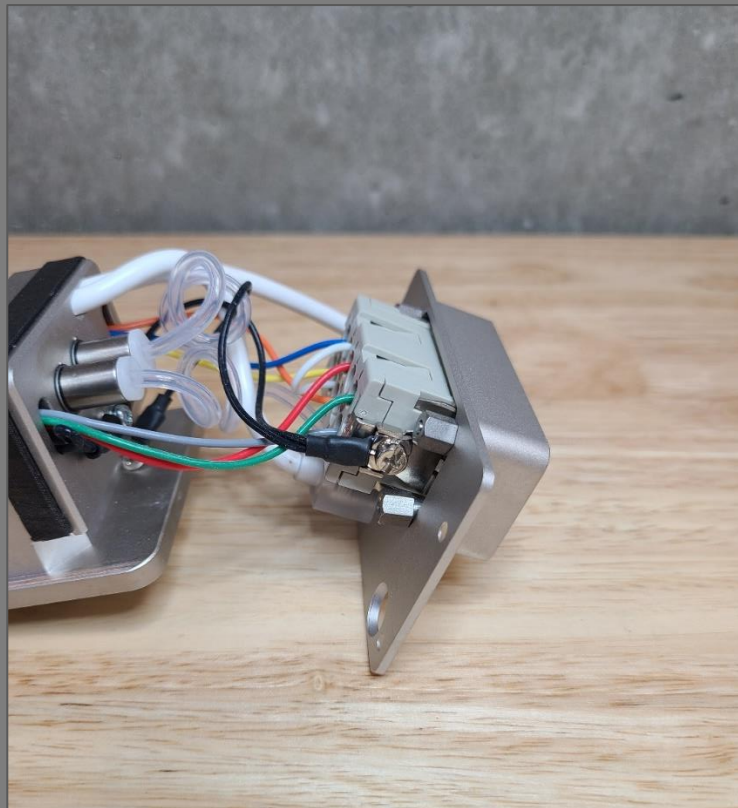

Separate the connector from the front plate.

2.3 Soldering Electronics

Add the nuts from the original Harting connector to the new one and mount the new connector to the front plate. Thread a 3/8-32 panel nut onto the HV plug to clamp it to the front plate and tighten it to prevent the HV plug from turning freely in the front plate.

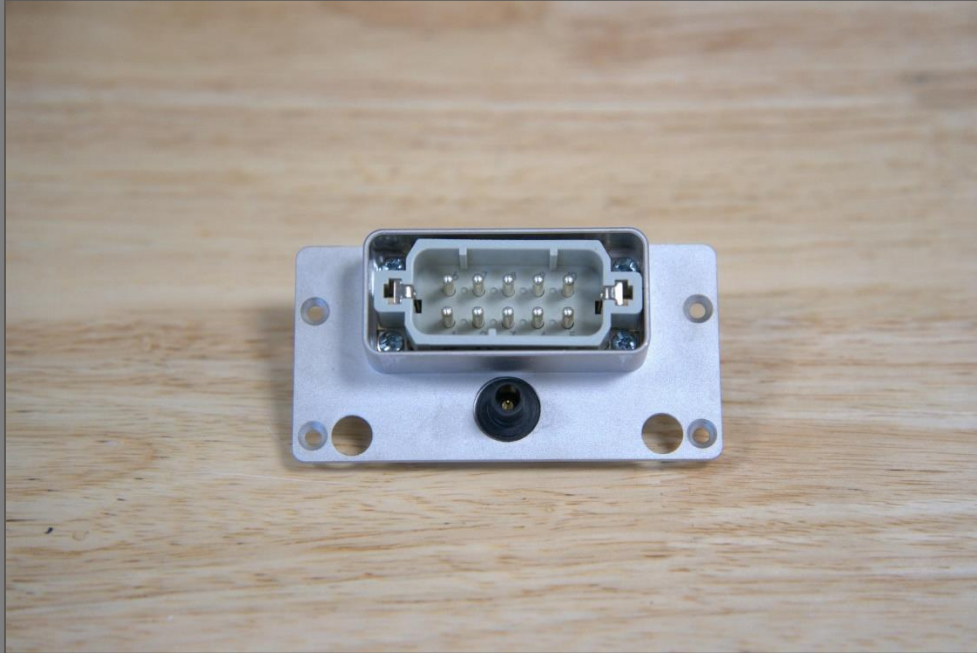


Cut a 10 cm section of HV wire and strip 1 cm from each end. Cut a second 60 cm section and strip 1 cm from one end. Solder the 22 MOhm resistor between the two sections of HV wire. Cut a 6 cm section of 1/4"- polypropylene tubing and slide it over the resistor.




Add heat shrink over the tubing. Add a short section of heat shrink to each end of the previous heat shrink to prevent the tube from sliding. Feed the long end of the HV cable through the port in the top of the 3D-printed electronics enclosure and carefully fold the resistor in.




Cut the plastic insulation off one end of a butt splice crimp connector using wire cutters.

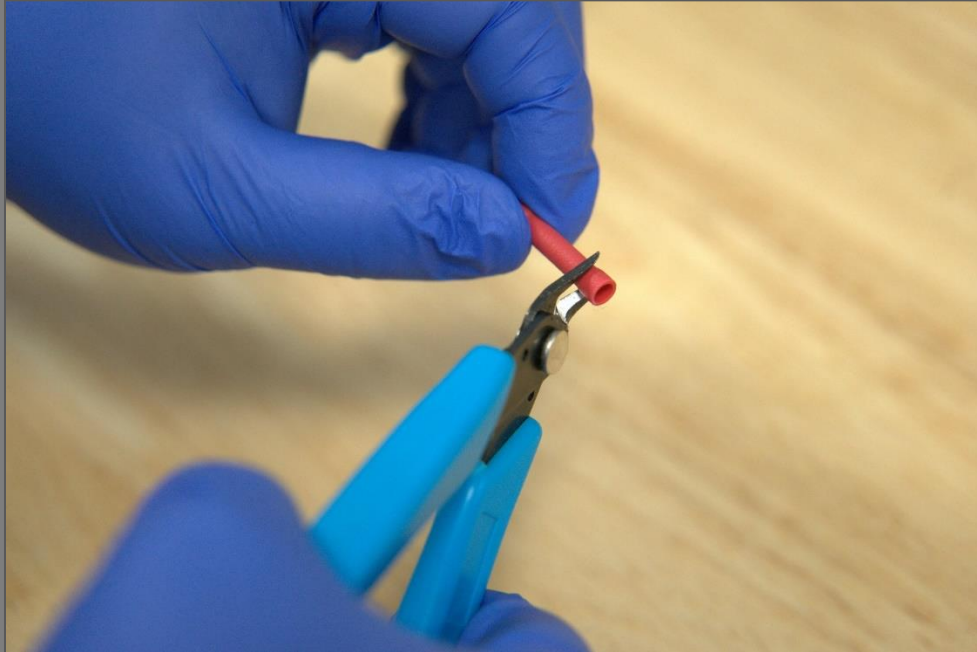

Slide the 3D-printed HV connector onto the short end of the HV wire with the threaded end facing out. Feed the short end of the HV wire into the uncut end of the crimp connector and crimp. Slide the cut end of the crimp connector over the gold-plated solder well on the inside end of the HV plug mounted to the front plate. It should slide on easily but make firm contact.

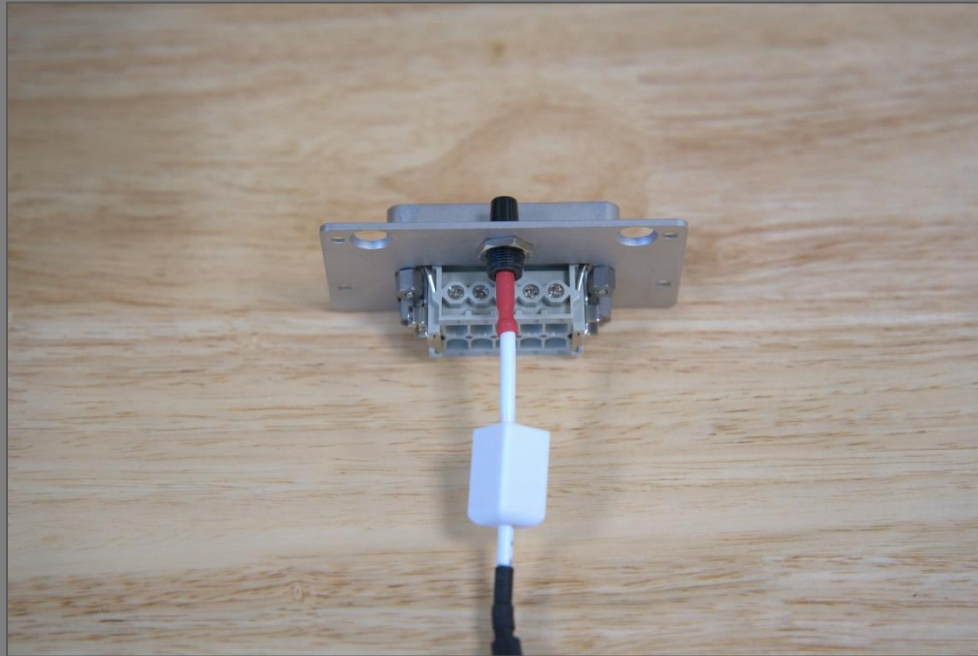

Heat shrink a 2.15 kOhm resistor, leaving the last 5 mm of each lead exposed. Strip 5 mm from each end of a 4 cm section of wire to make a jumper. Place the resistor between pins 6 and 10 in the Harting connector and place the wire jumper between pins 4 and 9. Tighten the set screws fix the connections.

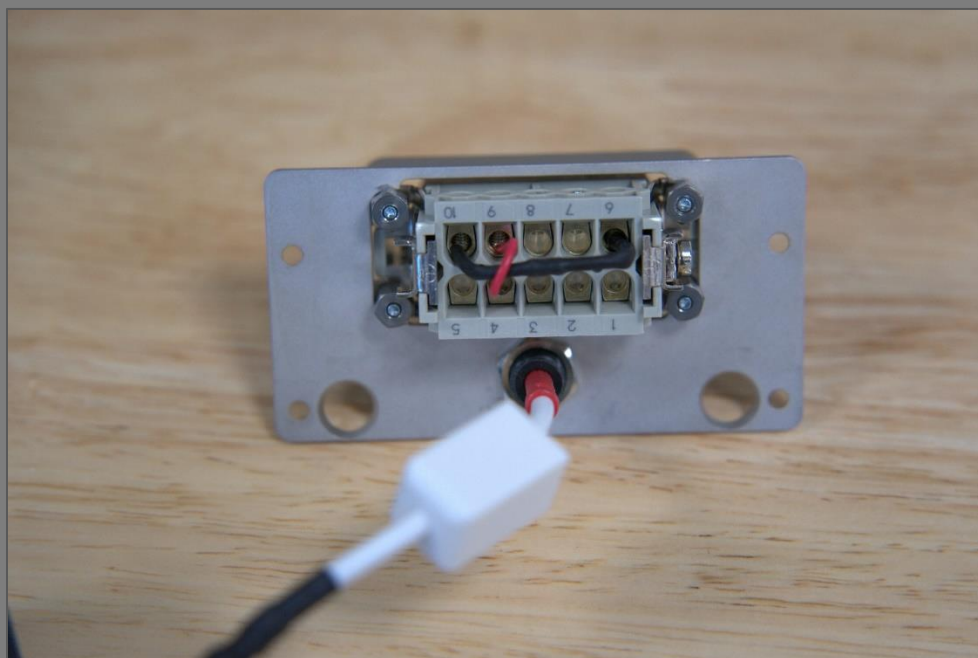

Thread the 3D-printed HV female connector onto the inside end of the HV plug and tighten to secure the crimped HV line to the HV plug. Carefully fold the HV cable into the 3D-printed electronics enclosure and mount the front plate to the electronics enclosure using the 2 of the original flush head mounting screws.



2.4 Assembling the shaft assembly


Cut the 1/4-40 lead screw to a length of approximately 9 cm using whatever cutting tool you have on hand. The threads on the cut end do not need to be preserved. Modify the lead screw nut so it fits into the shaft with the two mounting holes lining up with the corresponding threaded holes on the shaft. The nut is made of acetal/POM, so it is easily shaped with common tools (files, hacksaw, dremel, etc.)



Thread the uncut end of the lead screw into the lead screw nut from the bottom until approximately 29 mm protrudes from the top of the nut. Mount the lead screw nut in the 3D-printed shaft using two 6-32 screws.




Slide the shaft up through bottom of shaft guide by lining up the grooves in the shaft with the guide pins on the inside of the guide shaft. All 3D printers print slightly differently, so tolerances may need to be adjusted in the original design files. Variable sets have been included in each design file to simplify adjustments of critical dimensions. If the shaft does not slide smoothly in the shaft guide, you may need to increase the bore diameter of the shaft guide by editing the "Bore_diameter" variable in the FreeCAD file for the shaft guide. If the shaft is too loose in the shaft guide, decrease the bore diameter. If the shaft slides smoothly in the shaft guide but cannot pass the guide pins, you may need to increase the pin-to-pin distance ("Pin_to_pin_distance" variable in the FreeCAD file for the shaft guide) or decrease the pin radius ("Pin_radius" variable). If the shaft has rotational play, you may have to decrease the pin-to-pin distance or increase the pin radius.

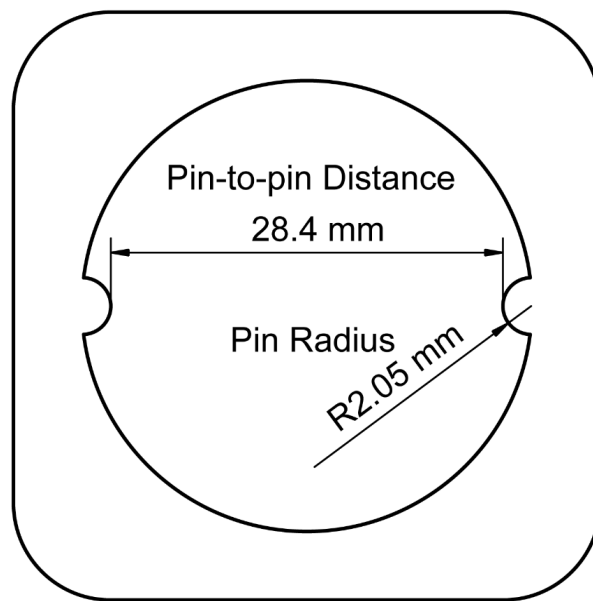

Push brass flanged sleeve bearing into the hole in the 3D-printed cap and press the M3 threaded insert into the side mounting hole from the inside.



Slide shaft collar over the top end of the lead screw and tighten the set screw 13 mm below the top.

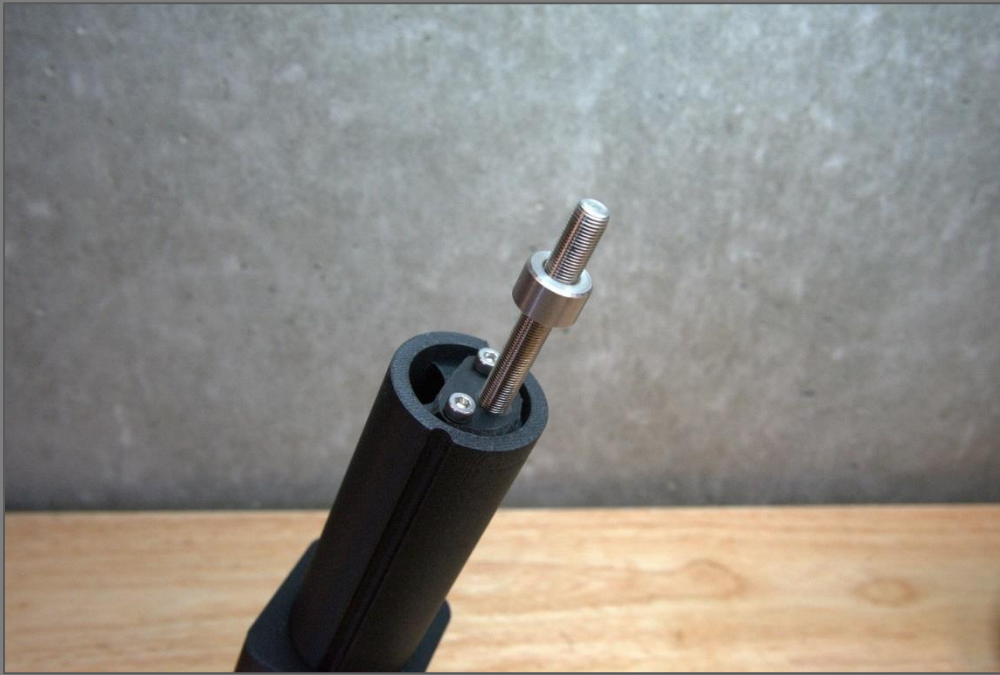

Guide lead screw up through the brass bushing in the cap. Slide the thumb knob onto the lead screw. While pushing the thumb knob against the shaft collar, tighten the set screw on the thumb knob to fix it to the lead screw. Turn the thumb knob clockwise while lining up the grooves on the shaft with the ridges on the inside of the cap until the cap slides over the shaft.



Thread M3 set screw into threaded insert on back side of cap, but don't tighten it all the way. Feed the free end of the HV wire coming from the electrical enclosure through the hole in the top of the cap and down through the shaft. Do the same with one end of the PEEK line.

2.5 Assembling the High-Voltage Termination.

Disassemble a Bic Round Stic® pen and flush the polypropylene ink cartridge tube out with methanol to remove all the ink. Cut a 60 mm section of the pen tubing. Strip 70 mm of insulation off terminus of HV wire. Slide a 60 mm section of heat shrink onto the stripped portion of the wire and heat.



Slide the plastic pen tube over the heat-shrunk portion of the wire all the way to the terminus of the insulation and solder the HV pogo pin onto the terminus of the wire.




Slide the plastic tube down onto the splines of the pogo pin and heat shrink where the plastic tube meets the HV wire insulation.




2.6 Installing the strip mount

(procedure is the same whether using the commercial strip cartridges or strip cartridges produced in-house).


Feed PEEK tubing and the HV terminus through the holes on the top surface of the cap and down through the shaft.

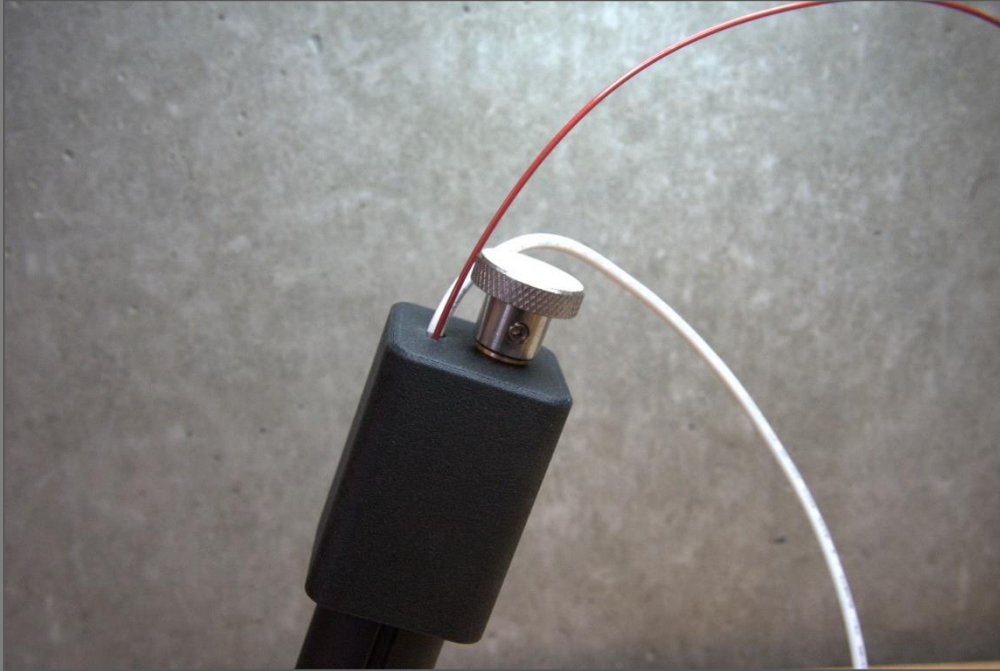

Slide a 10-32 threaded (male) 1/16-inch nut onto the end of the PEEK line. Slide a 10-32 to 1/4-28 adapter onto peak line and begin to thread the 1/16-inch nut into the adapter. Position the bottom of adapter 32 mm from terminus of PEEK line and tighten connection to clamp the PEEK fittings in place on the PEEK tubing.

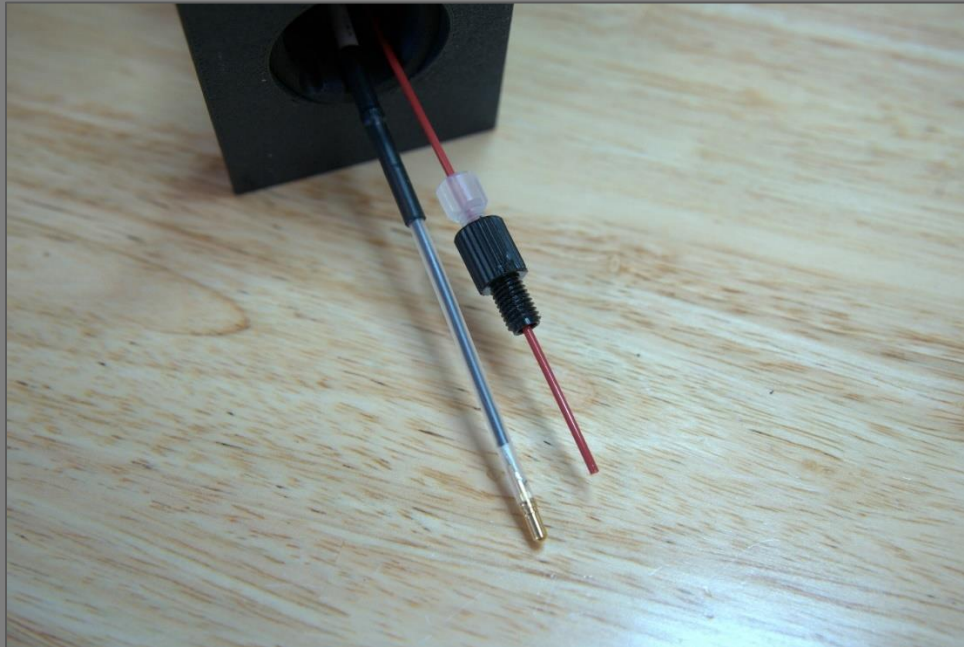


Feed the PEEK terminus into the threaded hole in the 3D-printed strip mount, and the HV terminus into the unthreaded hole. Thread in the 1/4-28 PEEK adapter while ensuring the PEEK tubing is rotating freely within the shaft assembly. Thread the 4-40 set screw into the threaded hole on the side of the strip mount but don't fully tighten.

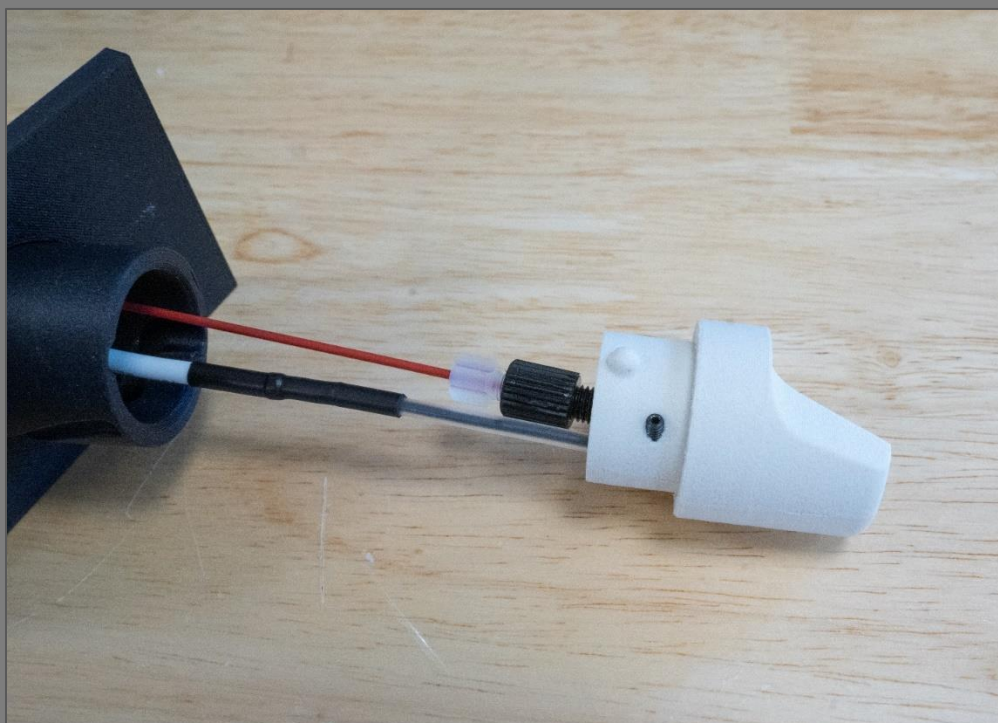

Insert a strip cartridge into the strip mount and ensure it is fully seated. If the strip cartridge does not slide in easily, you will need to increase the vertical tolerance by decreasing the “Strip_Height” variable in the FreeCAD file for the strip holder (if using opensource strip cartridges), or by increasing the “Strip_pocket_height” variable in the FreeCAD file for the strip mount (if using commercial strip cartridges). If the strip cartridge is loose in the strip mount, you will need to decrease tolerances. Slide the high voltage terminus into the strip mount until it contacts the rivet on the strip cartridge (look down along top of paper strip to ensure). Push it an additional mm to ensure good contact between pogo pin and rivet, but not so far that the pogo pin is fully depressed. Tighten set screw on the side of the strip mount to fix the height of the high voltage contact.



Turn the PEEK adapter nut clockwise until the PEEK terminus is roughly 1 mm above the surface of the paper strip (make sure the PEEK is turning freely within the shaft assembly so as not to introduce torsion in the PEEK tubing). The following image shows a strip mount with properly adjusted high voltage contact and solvent line.

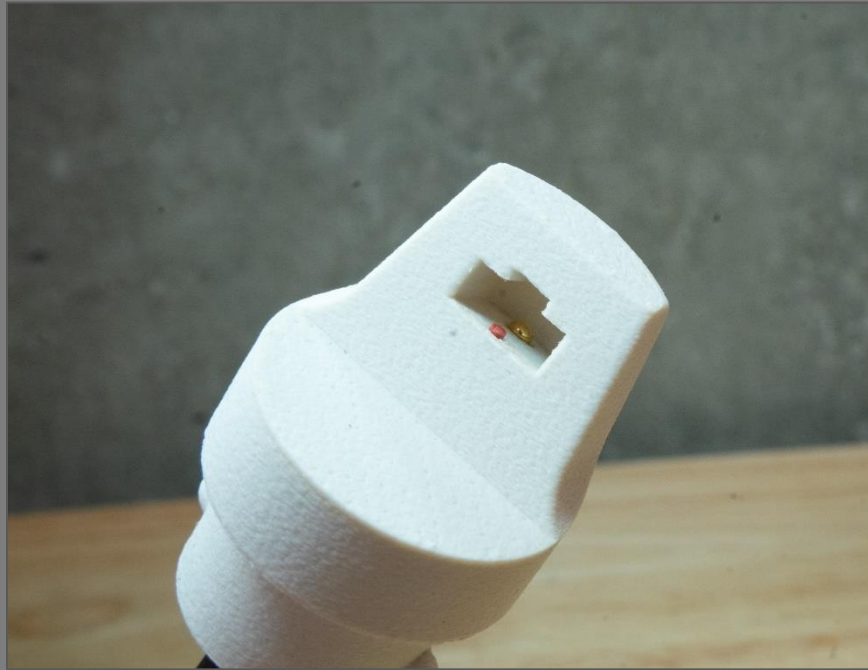

Line up locking lugs on insert with slots on the shaft, slide insert into shaft while twisting clockwise until insert clicks into place. Pull extra slack in high voltage line up through the cap.



2.7 Installing assembly in the ESI housing

Pass the electrical enclosure assembly and PEEK tubing up through the bottom of the positioning plates.



Slide the positioning plates over the shaft assembly and seat them on the shaft guide.



Mount the positioning plates and shaft assembly to the ESI housing and replace the 4 mounting screws. Replace the left thumbscrew only, to allow easier access to the loading port in the shaft guide when loading and unloading strip cartridges.




3.0 Opensource Strip Cartridge Construction

3.1 Cutting Paper Strips


If you choose to make your own paper strips, you will need to have a method of cutting the strips out of the chromatography paper. Paper strips can be cut by hand using a razor blade, by custom die cutters, or by a laser cutter. Whichever technique is chosen, it is crucial that the tip of the paper is cut to a fine point without any fraying otherwise you may experience unstable spray. If using a laser cutter, we recommend rinsing the strips in deionized water prior to use to remove pyrolysis products generated from the heat of the laser. We have included an .svg file for laser cutting. Laser settings must be optimized for your specific laser cutter.

3.2 Printing Strip Holders and Strip Holder Clips

The FreeCAD file for the strip holder contains the strip holder clip as well. Print the strip holders and the strip holder clips standing vertically with the back surfaces on the build plate. Support should not be necessary, as the overhangs are minimal.

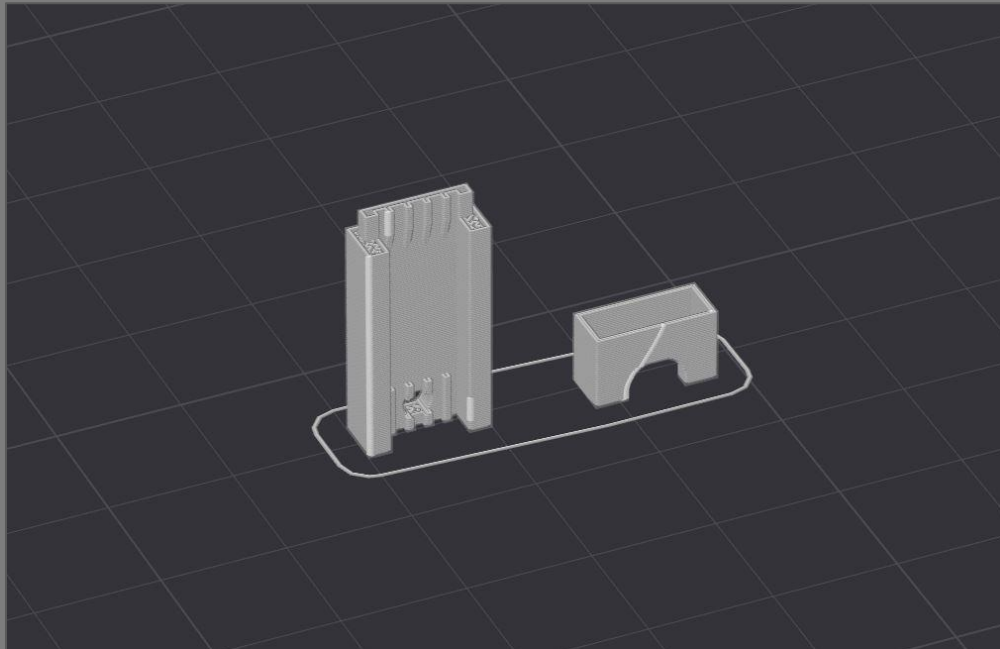

Print the strip tray laying flat. Support is recommended for the overhangs.

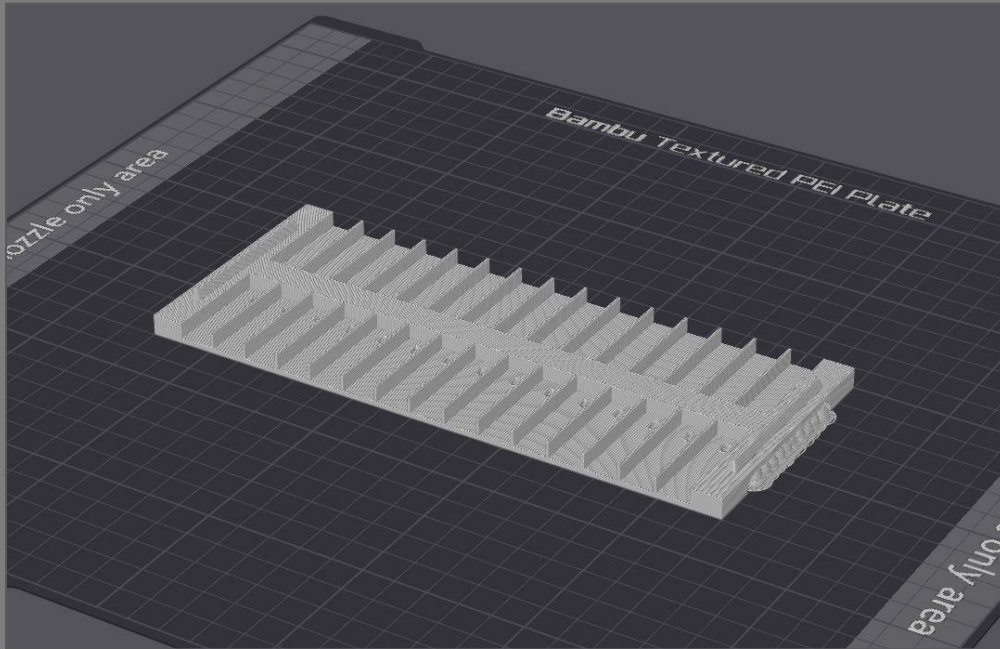


3.3 Assembling Strip Cartridges


Place a paper strip in a 3D-printed strip holder with the holes in the paper strip and strip holder aligned.



Insert an M2 rivet into the hole in the strip holder and press down until the paper strip is clamped flat against the strip holder.

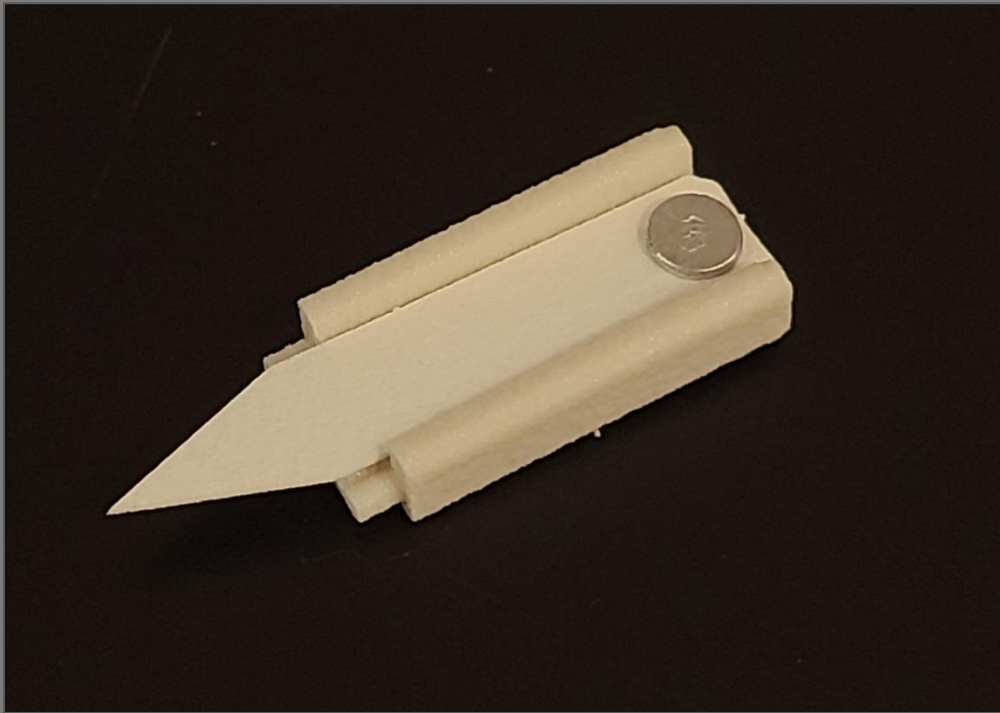

Slide a 3D-printed strip holder clip over the end of the paper strip and onto the strip holder until it clicks. Be careful not to damage the tip of the paper. The paper should be laying flat in the strip holder. You may need to modify tolerances in the FreeCAD file for the strip holder clip and/or strip holder to ensure a snug fit.

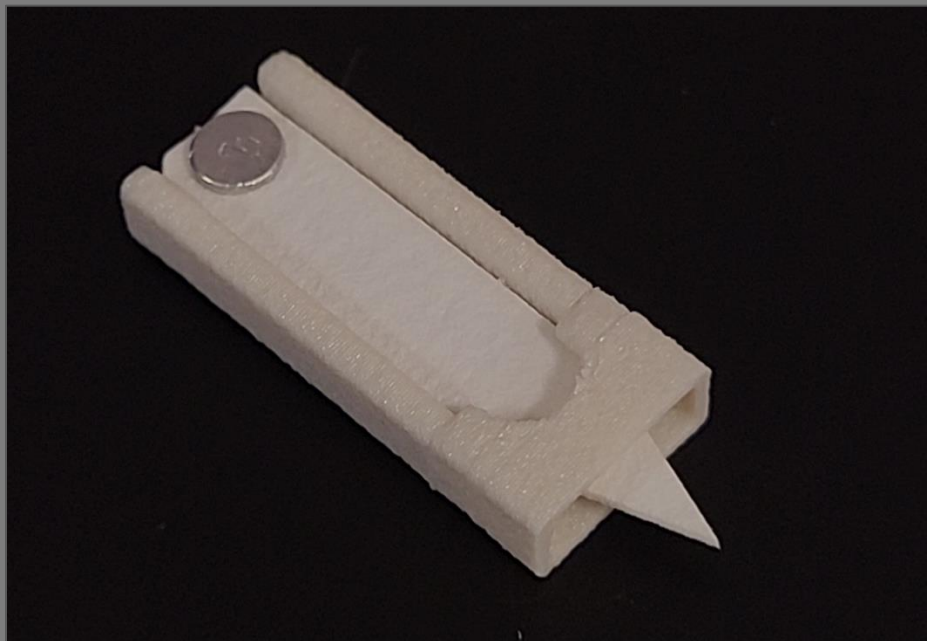


Put each assembled strip cartridge into a 3D-printed strip tray for storage until use.

4.0 Operating Instructions

The following instructions assume familiarity with instrument software and general instrument use. Refer to your instrument user manual for more specific instructions.

4.1 Positioning Calibration Procedure

- Turn off the heat to the ion transfer tube / inlet and allow it to fully cool.
- Loosen the 3 socket head screws and remove the left side panel from the source. Using your fingers, push the view port plug out of the source housing from the inside.

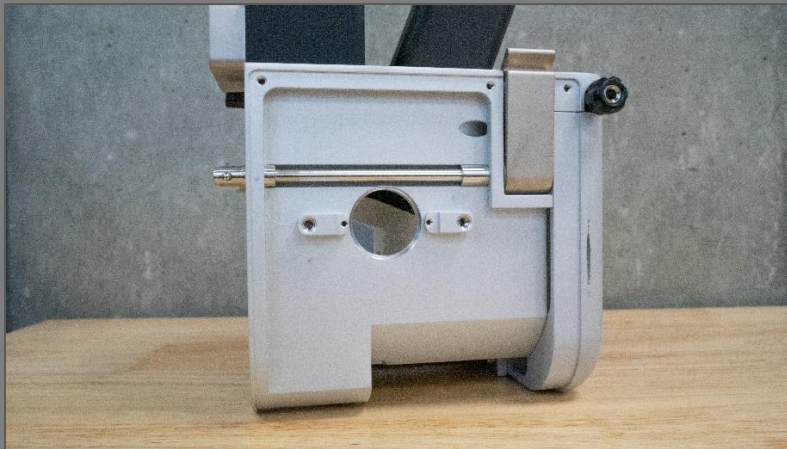

- With source removed from the instrument and the shaft at the bottom of its travel, insert a distance calibration strip and slide the shaft assembly to the back of its travel.
- Mount the source to the instrument.
- While looking through the side view port, slide the shaft assembly forward until the calibration strip is close to the inlet. Adjust the vertical positioning of the calibration strip by turning the thumb knob on the shaft assembly until the center of the calibration strip is aligned with the inlet. When the desired height is reached, tighten the set screw on the backside of the cap.



- To calibrate the gap between the inlet and paper strips, slide the shaft assembly forward until the calibration strip touches the inlet. It is crucial that the inlet is cool for this step, or the calibration strip may melt onto the ion transfer tube and contaminate the instrument. Tighten the black retainer knob on the source housing to clamp the shaft assembly in position.
- Remove the source from the instrument and remove the calibration strip from the strip mount. The viewport plug and side panel can be replaced.

4.2 Spray Solvent Delivery Set Up

The free end of the PEEK tubing can be terminated with a Luer-Lok adapter and connected to a 10 mL glass syringe filled with the desired spray solvent. We recommend priming the solvent delivery line before mounting the first strip cartridge to be run to remove air from the solvent line.

Using a syringe pump connected to the instrument allows solvent delivery to be programmed into the instrument method along with voltage and scan parameters.

4.3 Method Development

The source should be recognized by the instrument as an ESI source. Develop the instrument method as a direct infusion ESI method. If the source is not recognized, check the resistance across pins 6 and 10 on the Harting connector. You may need to replace the resistor if resistance is greater or less than 2.15 k Ω .

Spray solvent should be selected based on the matrix that your sample is in and the analytes of interest. Refer to the literature for examples of spray solvent composition for different applications.

Program solvent delivery to start at the beginning of the method before voltage turns on to give the analytes time to dissolve into the solvent and the solvent time to wick to the tip of the paper. It takes approximately 120 μ L of spray solvent to wet a paper strip, so the solvent flow rate will determine the amount of time needed to pre-wet the paper. Exact solvent flow rates and deposition volumes can be altered to optimize sensitivity and achieve stable spray for the desired method duration. The amount of solvent on the paper can affect the spray mode, which in turn affects ionization and transmission efficiency. A good starting point is a flow rate of 300 μ L/min for 0.4 minutes

Program the voltage to come on once the paper strip is wet. We recommend optimizing spray voltage for the spray solvent, sample matrix, and analytes you are interested in, but 3600 – 4000 V are typical for paper spray methods in positive mode. Negative-mode paper spray requires lower spray voltages to prevent corona discharge and charring of the paper. The duration of high voltage depends on the desired length of the method. Paper spray methods can range from a few seconds to several minutes.

Program data acquisition to continue for several seconds after the voltage turns off so the signal returns to baseline at the end of the method. This will simplify peak integration during data processing.

4.4 Sample Preparation

A diverse variety of sample types are amenable to paper spray mass spectrometry including but not limited to biofluids, water samples, health supplements, pharmaceuticals, forensic residues, and tissue homogenates. In most cases, little to no sample preparation is required and samples can be directly deposited on the paper strips without prior filtration or extraction steps. A sample deposition volume of 10 μL is standard, but lower volumes are sometimes used for complex sample matrices.

Samples should be allowed to dry on the paper, either at ambient conditions or with the help of a drying oven or desiccator prior to analysis.


4.5 Loading and Analyzing a Strip Cartridge

To load a strip cartridge, raise the shaft to its upward position by pulling upwards on the cap. The shaft and cap rotate 90° clockwise in the shaft guide to allow strip cartridges to be loaded from the cutout in the right side of the shaft guide.



Ensure the strip cartridge is properly seated in the strip mount and slide the shaft to the bottom of its travel. The bottom of the cap should be flush with the top of the shaft guide.

The strip cartridge is now loaded, and you can start the instrument method. You should get a sharp increase in signal when the voltage turns on and a sharp decrease in signal when the voltage turns off.

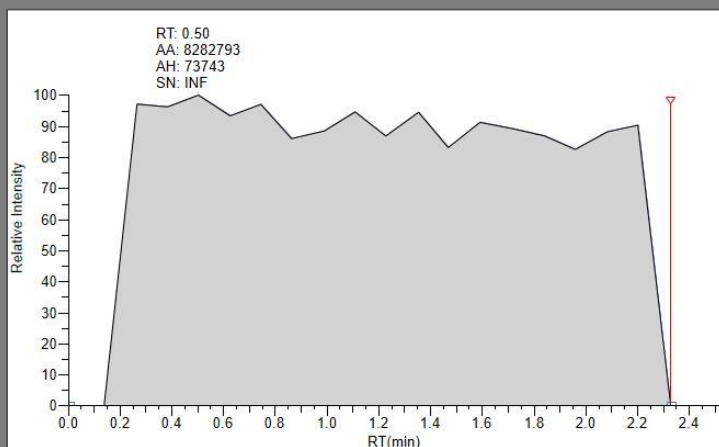

4.6 Troubleshooting

A common issue is solvent overflow caused by excess solvent deposition. This will be noticeable when changing strip cartridges as the sides and bottom of strip holders will be wet with solvent. When this occurs, reduce either the solvent flow rate or the duration of deposition.

If the signal doesn't come on, there is either a problem with the instrument method, the electrical contact, the solvent delivery, or the positioning of the paper strip relative to the inlet. Try the following steps in order:

With the voltage turned off, raise the shaft to the upward position and visually inspect the paper strip for solvent. It should be saturated with solvent but not overflowing or forming a drop at the tip. If there is no solvent, prime your syringe pump to remove air and check your instrument method to confirm it is programmed correctly to deliver solvent. You may need to increase the solvent delivery rate / duration.

If solvent is flowing but the paper is not getting wet, solvent may be wicking up the outside of the PEEK tubing into the strip mount. To resolve this, increase the inner diameter of the PEEK tubing.

If the paper looks saturated with solvent and there is still no signal, try moving the shaft assembly forward to decrease the inlet gap.

If there is still no signal, remove the source from the instrument and use a multimeter to confirm continuity between the pogo pin and the HV plug. There should be approximately 22 M Ω of resistance. If there is no continuity, remove the electronics enclosure and determine where the connection is broken.

If the HV line has appropriate resistance, check your instrument method to confirm the voltage and scan parameters are correctly programmed.

If the signal onset is delayed or choppy, try increasing the voltage onset delay or move the shaft assembly forward or back.

4.7 Unloading Used Strip Cartridges

To unload a strip cartridge after data has been acquired and the high voltage is turned off, raise the shaft to the upward position and use gloved fingers or forceps to extract the used strip cartridge. If there is solvent in the strip well on the strip mount, lower the solvent deposition volume for subsequent runs to prevent carryover. The high-voltage contact rivets can be washed and re-used, but the rest of the strip cartridge should be properly disposed of.